

General Information

Intended use:

Pancreatic Elastase Turbilatex® Combo is a latex turbidimetric assay **only for the quantitative detection of Pancreatic Elastase E1 in human stool samples** (not to be used for body fluid as blood, serum, plasma, urine, cerebrospinal fluid, oral fluid, synovial fluid or empyema fluid). This assay is simple and widely applicable. This product is optimized for several automated analyser.

For professional *in vitro* diagnostic use only.

Reagents:

Materials provided by CerTest Biotec:

Reagents	Code
Turbidimetric reagents (R1 & R2)	TL-022EL01 TL-022EL02
Auxiliary Reagents	
Calibration kit	TL-022EL70, TL-022EL71 TL-022EL72 TL-022EL73 TL-022EL74 TL-022EL75
Controls kit	TL-022EL08 TL-022EL09
Sample dilutions vials	MST-0019U

Preparation of reagents: .

R1 and R2 are ready to use.

Calibrators are ready to use.

Controls are ready to use

Storage and stability:

Kit components must be stored at temperature indicated on the label. Do not freeze.

Reagents are stable up to the expiration date printed on the label, always considering that reagent containers must be properly closed to avoid any contamination, must be kept away from the sunlight and conserved at temperature indicated on the label of each reagent.

Specimen:

Collect enough quantity of human solid stool samples. These samples should be collected in clean and dry containers (no preservatives or transport media). The samples can be stored in the refrigerator (2-8°C) for 7 days prior to testing. If not immediately tested, freeze the stored samples at -20 °C maximum 6 months. In this case, the sample will be totally thawed and brought to room temperature (15-30°C) before testing. Homogenize stool samples as thoroughly as possible prior to preparation.

Use **Pancreatic Elastase Turbilatex® Combo** stool collection tubes for sample collections described the instructions for use.

Assay procedure

Pancreatic Elastase Turbilatex® Combo can be performed on every open chemistry analyser. Please follow the subsequent instructions in order to assure performance characteristics as describes in the instructions for use. This instruction has been validated by CerTest BIOTEC S.L Laboratories.

Additionally, please read the “Instructions for use” for instructions on operating and programming user defined test.

Application parameter set up:

Specific analyzers settings for **Pancreatic Elastase Turbilatex® Combo** must be programmed onto the analyzer, see below. For instructions, consult the TC220 (Tecom) analyzer manual and instructions for use provided with the kit.

Loading of reagents:

Load reagents according to the TC220 (Tecom) analyzer manual.

Calibration curve establishment:

A 6-points calibration curve can be established in TC220 (Tecom) analyzer. For instructions consult analyzer manual.

Calibration stability:

Calibrate the system at least once a month is extremely recommended. Recalibrate the system when reagent lot is change or when the controls are out of the assigned range given in the control label and CoA.

QC controls:

Pancreatic Elastase Turbilatex® Combo controls C1 and C2 must be assayed each day before running patient fecal sample extract to validate the calibration curve. The controls have assigned value ranges indicated on the label and certificate of analysis supplied. The control measurements must be within the indicated value range to obtain valid results for patient fecal extract. If the control values are out of range, follow next procedures: 1) Repeat QC control measurement, 2) Repeat calibration measurement.

Results:


The results are evaluated automatically by the analyzer and presented in µg hEL/g of stool.

Performance characteristics (*)

The following results have been obtained during the validation of **Pancreatic Elastase Turbilatex® Combo** on the TC220 (Tecom) analyzer.

Linearity:

Pancreatic Elastase Turbilatex® Combo using calibrator kit is linear in the calibration range of 18.8-400 µg hEL/g of stool.

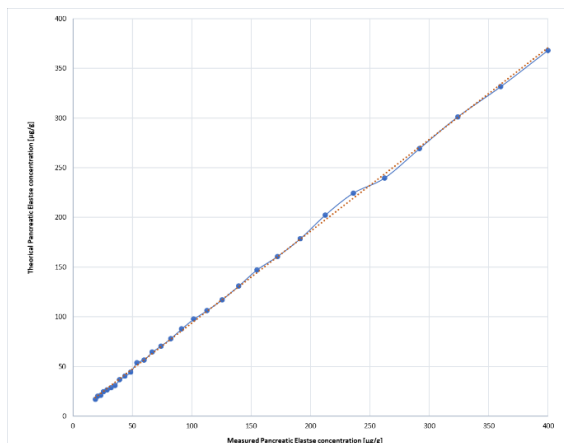

Measuring range:

Pancreatic Elastase Turbilatex® Combo assay measuring range is 5.3-625 µg hEL/g of stool. Samples more concentrated than 400 µg hEL/g of stool must be diluted for proper quantification by the user, using additional sample buffer.

Prozone effect:

Studies have been made up to a concentration of 10000 µg hEL/g of stool and no false negative results have been observed. Studies using higher concentrations have not been carried out. Samples with concentrations up to 625 µg hEL/g of stool can be measured without inhibitory prozone effect.

Detection limit:

Limit of detection (LOD): 1.07 µg hEL/g of stool. The lower limit of detection of **Pancreatic Elastase Turbilatex® Combo** was determined on 20 samples and 2 sample replicates as the mean value +2· SD.

Limit of quantification (LOQ): 5.3 µg hEL/g of stool. The lower limit of quantification is defined as the lowest actual amount of analysis that can be reliably detected; imprecision is < 20% as CV%.

Precision:

Pancreatic Elastase Turbilatex was tested with three different controls levels.

	Low (25 µg/g)	Medium (100 µg/g)	High (400 µg/g)
N	20	20	20
Mean (µg/g)	25.0	101.9	390.6
SD (µg/g)	0.8	1.7	9.4
CV (%)	3	2	2

Method comparison

The following tables provide information on the summary of the evaluations that have been carried out on the product **Pancreatic Elastase Turbilatex® Combo** with a cut-off of 200 µg of hEL/g of stool:

Pancreatic Elastase Turbilatex® vs Evaluation criteria		
	Mean Value	95% confidence interval
Sensitivity	89.1%	77.8-95.9%
Specificity	98.4%	94.4-99.8%
PPV	96.1%	86.5-99.5%
NPV	95.4%	90.2-98.3%
LR+	56.13	14.15-222.7
LR-	0.111	0.052-0.236

Shipping damage

Please notify your distributor, if this product was received damaged.

Symbols key

	For in vitro diagnostic use only		Keep dry
	Consult instructions for use		Temperature limitation
	Catalogue number		Lot number
	Use by		Manufacturer
	Contains sufficient for <n> test	DIL	Sample diluent
	Keep out of the sunlight		

Manufacturer

CERTEST BIOTEC

Pol. Industrial Río Gállego II, Calle J, N° 1, 50840,
San Mateo de Gállego, Zaragoza (SPAIN)
www.certest.es

NOTES

Please refer to the instructions for use for the detailed information about the test on the following:

Synthesis; Principle; Precautions; Reagents; Specimen collection; Interpretation of results.

(*) Analytical performance data were obtained with the Biolis 24 i(Tokyo Boeki) analyser.

TC220, Tecom/ Application parameters

ASSAY PARAMETERS	
Std. No	6
R1	150 µL
Sample	10 µL
R2	15 µL
Others	N/A
Reaction mode	Two points
Primary wavelength	510 nm
Secondary wavelength	None
Direction	Increase
Reagent blank lecture (cycle)	10 s after R2 addition
Final lecture (cycle)	300 s after R2 addition
Reaction time	close to 10 min
Linear range	18.8-400 µg/g
CALIBRATION	
Calibration Method	Spline
Calibration set	5 calibrators + Blank
Blank	Calibrator 0 (0 µg/g)
Calibrator 1	Calibrator 1 (25 µg/g)
Calibrator 2	Calibrator 2 (50 µg/g)
Calibrator 3	Calibrator 3 (100 µg/g)
Calibrator 4	Calibrator 4 (200 µg/g)
Calibrator 5	Calibrator 5 (400 µg/g)
STEPS	
Addition R1	0 s
Addition Sample	120 s
Incubation R1+S	120-180 s
Addition R2	220 s
Blank Lecture	230 s
Incubation reaction	300 s
Final lecture	530 s

TC Series Automatic Chemistry Analyzer----- Temp. Cubetas:36,94

Test Informes QC AJUSTES Detener Salir

Test Bioq. Otros Parametros Param. Calculados Ajuste Perfiles Secuencia Ensayo Secuencia Impresion

Item: CALP2, TF, HPylori, FOBh, **ELASTASA**

Simbolo: ELASTASA Nombre: Elastasa pancreatica Cod. de barras Param.Basicos Param.Calibracion

metodo: Two Points 1ª Long.O: 405 2ª long.O: NO Decimal: 1 Und: ug/g (Rtvo.1+Rtvo.2)+S

Rtvo.1 Vol.(10-400ul): 150 Incub.Rtvo: 60 Rtvo.1 Pos: 6 Vol. Del envase: 20 ml 2ª posicion: 0

Rtvo.2 Vol.(0-400ul): 100 Rtvo.2 Pos: 7 Vol. Del envase: 20 ml 2ª Posicion: 0

S Vol.(1.0-50.0ul): 3,0 Incub.Rtvo: 20 Vel. de Mezcla: 0 Tiempos de reaccion: 10 Check Time(s): 360

Rango Linealidad: 0 — 400 Correccion de resultados: 1,0000 X + 0,0000 ratio de Dilucion (2-150): 10

Absorbance Range: -5,0 — 10,0 Agotamiento de Sustrato: < 10,0 Vol. Muestra (3-70): 15

Lim. de Linealidad: 10 de Juicio de valor Alto Discrepancia: 10 ctado en la Pendiente = 0 Vol. Total (120-500): 150

Blank read Point offset: 0 Ratio de Dilucion: 1:9

test blanco Reactivo Blanco de Rtvo.: 0,0 — 10,0 Blanco Reactivo: 0,3765

Fabricante: Certest Biotec No. Lote: 200EDELRU00 Caducidad: 31/01/2021

Añadir Borrar Guardar

Super User Nombre:aaa 24/02/2020 15:16:24 Vers

Informes QC AJUSTES Detener Salir

Otros Parametros Param. Calculados Ajuste Perfiles Secuencia Ensayo Secuencia Impresion

Simbolo: ELASTASA Nombre: Elastasa pancreatica Cod. de barras Param.Basicos Param.Calibracion

Nº Std: 6 Nº de Rept.: 2

Calibracion: Spline

Lot: 67EL004

Caducidad: 31/01/2021

Dilucion

K: 0,0 El orden de posicion de los Stds. Deberá ser de MENOR a MAYOR, de lo contrario no se calibrara

R0: 0,000872728

a: 0,000552106

b: -0,000000954

c: -0,000000003

d: 0,0 Hª de Calib Recalc.

	1	2	3	4	5	6
Posicion Std	1	2	3	4	5	6
Valor Std	0	25,0	50,0	100,0	200,0	400,0
Absorbancia	0,0009	0,0140	0,0254	0,0406	0,0671	0,0933

Añadir Borrar Guardar

Super User Nombre:aaa 24/02/2020 14:11:10 Vers